




Nulife
I.V. Cannula

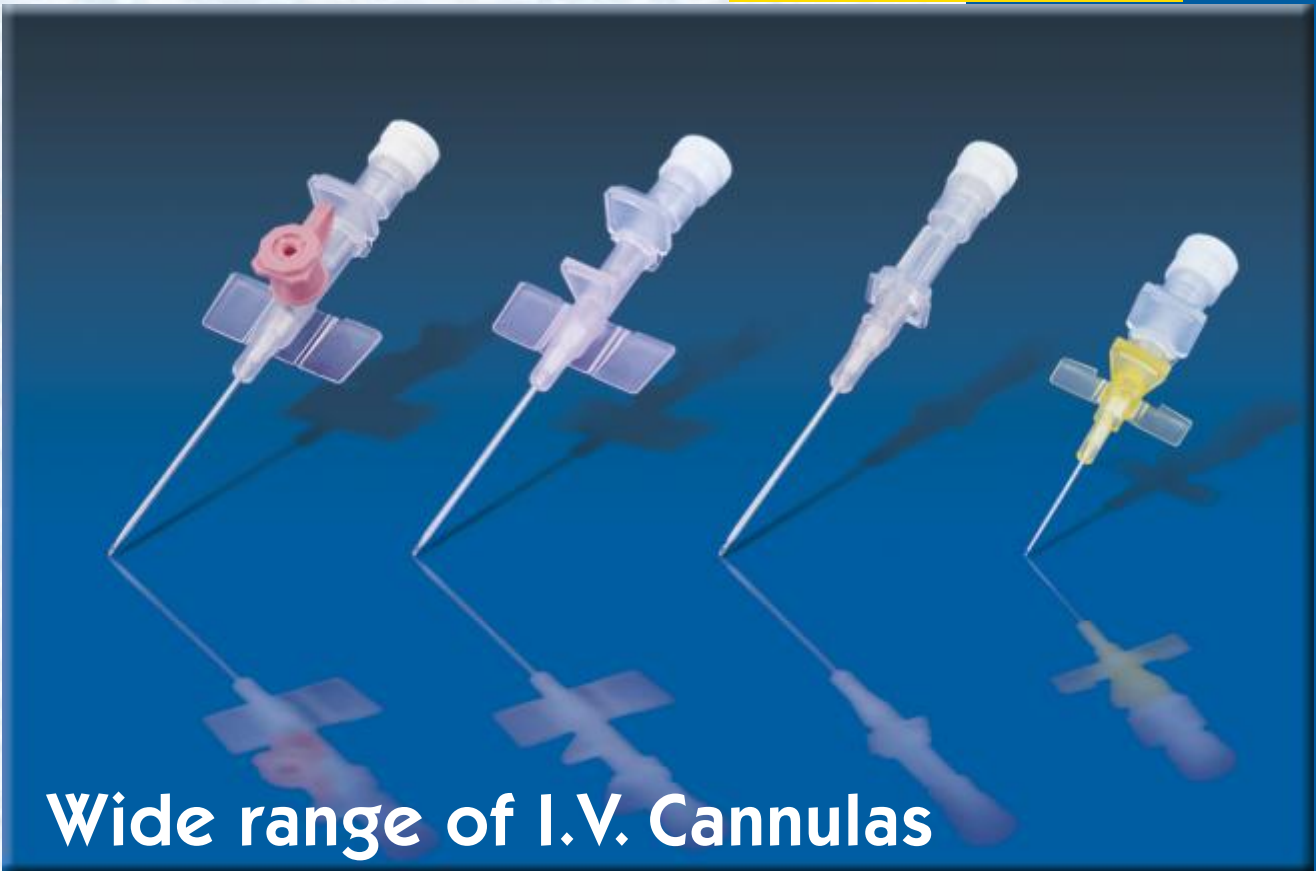


Wide range of I.V. Cannulas

for every use.



NULIFE and **YOU**
Together let's save lives

Nulife I.V. Cannula

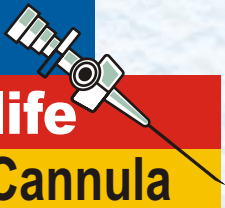

USER FRIENDLY I.V.CANNULAS FROM NULIFE

Special Features:

- Specially tapered **ptfe** catheter is a precision formed to perfection and **kink** resistant for easy vein puncture with minimum trauma.
- Cannulas have a Japanese needle and luer lock Plug.
- Sterile, disposable and non-pyrogenic.
- The angled flexible wings allow for easy fixation which prevents pistoning or rolling of cannula on the patients body.
- Latest tipping technology with special emphasis on the Teflon tip for easy vein penetration with minimal patient trauma.

NUCATH I.V. Cannula for vein cannulation in pediatric and geriatric patients

- Miniature design to deal with specific problems.
- Small projection on main body.
- Small flexible wings.
- No injection valve or side port.

NUFLON I.V.Cannula with injection valve and wings

- Injection site with one-way retracting valve.
- Allows needle-free injection of medicine, reducing chances of infection and needle stuck injuries.
- Specially designed valve cap incorporating a recessed plug with protective skirt to effectively prevent contamination when valve is not in

Also Available:

- **NUMED M I.V. Cannula** without port & without wings.
- **NUMED R I.V. Cannula** with small wings without injection valve.
- **NUMED K I.V. Cannula** with injection valve & suturable.



Gauge rate	Colour	Catheter O.D.	Catheter length	Water flow
			in mm	ml. /minute
14	Orange	2.1	45	300
16	Grey	1.7	45	200
17	White	1.5	45	140
18	Green	1.3	45	90
20	Pink	1.1	33	61
22	Blue	0.9	25	36

Packaging: Every I.V. Cannula is individually packed in rigid blister pack. 50 packed I.V. Cannula are available in a dispenser box.



MRK HEALTHCARE

B4/B5, Byculla Service Industries, D. K. Marg,
Mumbai 400027. Tel: 91-22-23748371/2/3.
Fax: 91-22-23748370.
E-mail: mrkhealthcare@vsnl.com
Website: www.mrkhealthcare.com

